



Sentri™ end tidal CO₂ monitoring range

Masks and nasal cannula



Oxygen and Aerosol Therapy ▪ Variable Oxygen Concentration (Low Flow)



Sentri™ end tidal CO₂ monitoring range

Mask or nasal cannula the choice is yours

The Sentri product range is available as an adult mask and in three sizes of nasal cannula. Both options permit the sampling of exhaled carbon dioxide in non-intubated patients during the administration of supplementary oxygen.

There are also options available for connection to Microstream® capnography.

[Find out more](#)

Nasal cannula available in three sizes

To accommodate a wide range of patients from infant to adult

Soft, curved prongs

For improved anatomical fit and comfort

Oxygen delivery and ETCO₂ sampling

Non-slip toggle

Provides a secure fit



Oxygen is delivered through one nasal prong and exhaled gas is sampled from the other. Nasal cannula may be more appropriate for paediatric patients when high oxygen flows may “dilute” the CO₂ sample and give a low (or no) value. A face mask may be more appropriate when the nares are occluded or obstructed.



Sentri™ Intersurgical EcoLite™ mask kit

Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, whilst the second softer TPE material is utilised in the seal, which is in contact with the patients face.

Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilisation of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask body and seal, reducing the environmental impact.

Intersurgical EcoLite

[Find out more](#)

Eco range

[Find out more](#)



Incurved nose seal

Conforms to different nose shapes designed to minimise the risk of oxygen entering patient's eyes

Multi-channel oxygen tube

Oxygen still flows even if tube is kinked

Luer lock port

Secure ETCO₂ monitoring line connection

Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

Soft face seals

Increased patient comfort




A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

Two integral chin seals ensure the mask fits a wider range of patient face shapes



Products available with CO₂ monitoring line for Microstream® capnography

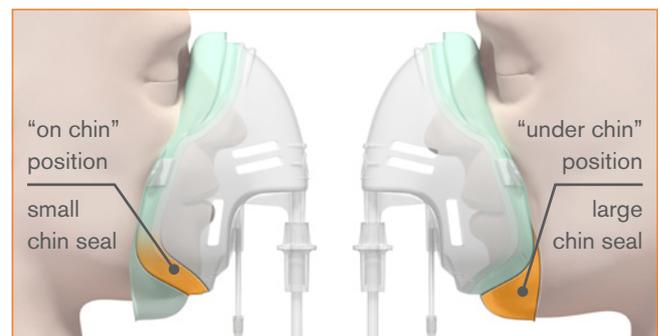



Available with oxygen tubing featuring a universal connector, which can be attached to a 6mm flow meter spigot or fitted directly onto a 9/16-18 UNF Diameter Index Safety System (D.I.S.S.) flow meter



Senti™ Intersurgical EcoLite™ mask

| Code | Description | Tube length | Box Qty. |
|---------|---|-------------|----------|
| 1142015 | Senti Intersurgical EcoLite, adult, mask with CO ₂ monitoring line and tube | 2.1m | 30 |
| 1143015 | Senti Intersurgical EcoLite, adult, mask with CO ₂ monitoring line, filter and tube | 2.1m | 30 |
| 1143018 | Senti Intersurgical EcoLite, adult, mask with CO ₂ monitoring line, filter and tube with universal connector | 2.1m | 30 |
| 1141015 | Senti Intersurgical EcoLite, adult, mask with tube | 2.1m | 30 |
| 1142016 | Senti Intersurgical EcoLite, adult, mask with CO ₂ monitoring line for Microstream® capnography and tube | 2.1m | 30 |
| 1143017 | Senti Intersurgical EcoLite, adult, mask with CO ₂ monitoring line for Microstream® capnography, filter and tube | 2.1m | 30 |

Make an enquiry



Senti nasal cannula

| Code | Description | Tube length | Box Qty. |
|---------|--|-------------|----------|
| 1144001 | Senti, adult, nasal cannula with curved prongs and tube | 2.1m | 50 |
| 1144002 | Senti, adult, nasal cannula with curved prongs, CO ₂ monitoring line, filter and tube | 2.1m | 40 |
| 1144005 | Senti, paediatric, nasal cannula with curved prongs and tube | 2.1m | 50 |
| 1144006 | Senti, paediatric, nasal cannula with curved prongs, CO ₂ monitoring line, filter and tube | 2.1m | 40 |
| 1144009 | Senti, infant, nasal cannula with curved prongs and tube | 2.1m | 50 |
| 1144010 | Senti, infant, nasal cannula with curved prongs, CO ₂ monitoring line, filter and tube | 2.1m | 40 |
| 1144015 | Senti, paediatric, nasal cannula with curved prongs, CO ₂ monitoring line for Microstream® capnography, filter and tube | 2.1m | 40 |
| 1144016 | Senti, infant, nasal cannula with curved prongs, CO ₂ monitoring line for Microstream® capnography, filter and tube | 2.1m | 40 |
| 1144017 | Senti, adult, nasal cannula with curved prongs, CO ₂ monitoring line for Microstream® capnography, filter and tube | 2.1m | 40 |

Make an enquiry

Illustrated Lower environmental impact product

Microstream is a registered trademarks of Oridion Medical 1987 Ltd.

IS10.24 • Issue 16 03.24



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 965 6300 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 14001:2015, ISO 9001:2015, ISO 13485:2016 and MDSAP

Please think before you print
Save energy and paper.
If you must print this information sheet please print it double sided.